TECHNICAL DATA

V2 VAULTZ

V03 | 13.08.2026

Table of Contents

1. Base components	4
2. Technical specifications	4
2.1 Cascade refrigeration unit, air-cooled; alternatively water-cooled	4
2.2 Electrical control cabinet, machine	5
2.3 Electrical control cabinet, cell	6
2.4 Special evaporator in the anteroom	6
2.5 Special evaporator in the treatment room	6
2.6 Cell combination TYPE ISO 150/160	7
2.6.1 All-glass entrance door	8
2.6.2 Intermediate door	8
2.6.3 Deluxe window to the treatment room	8
2.7 Service router with interface for remote monitoring	9
2.8 Camera in the treatment room	9
2.9 Music playback in the treatment room	9
2.10 Communication system from control unit to treatment room	9
2.11 Operating unit and additional operating options	9
2.12 Handrail in the treatment room	10
3. Optional components	10
3.1 Water-cooled condenser	10
3.2 Air-cooled condenser	11
4. Optional components	11
4.1 Hybrid cooling	12
4.2 Emergency Door	12
4.3 Wood panelling incl. handrail in the treatment room	12
4.4 Chiller	12
4.5 Special dimensions V2 cell	13
4.6 Big condenser	13
4.7 Camera in the anteroom	13
4.8 Music playback in the anteroom	13
4.9 Heat recovery	14
4.10 Division of the cooling unit	14

4.11	Frost heave protection	14
5.	On-site services	15
5.1	Connections	15
5.2	Set-up	15
5.3	Ambient conditions	16
6.	General provisions	18
7.	Pictures	19
8.	Certificates	20

TECHNICAL DATA V2 VAULTZ

These technical specifications are original and have been prepared by the manufacturer.

Any reproduction requires the consent of: L&R Kältetechnik GmbH & Co. KG

Sundern, 13.08.2026, Version 03

1. Base components

- Machine part, including machine housing, designed as a cascade refrigeration unit, air-cooled (alternatively water-cooled)
- Electrical control cabinet, machine
- Electrical control cabinet, cell
- Special evaporator in the anteroom
- Special evaporator in the treatment room
- Cell combination TYPE ISO 160 in its two-room version
- Service router with interface for remote monitoring
- Camera in the treatment room
- Music playback in the treatment room
- Intercom system in the treatment room
- Control panel and optional operating options
- Handrail in the treatment room

2. Technical specifications

2.1 Cascade refrigeration unit, air-cooled; alternatively water-cooled

- Highly efficient machine design
- Special heat-resistant and pressure-resistant heat exchanger plates as cascade intermediate coolers
- Additional pre-heaters to reduce the temperature of hot gas
- Control of all stages with electronic expansion valves
- Highly efficient oil management system
- All necessary and relevant safety components for monitoring and shutting down individual components for maintenance-friendly operation of the system

Cooling capacity approx.:	-6.9 kW
umber of refrigeration circuits:	3 units
Refrigerants:	
1. Stage for the condensation of the 2nd stage of cascade	R 449 A
2. Stage for cooling the anteroom and condensation of the 3rd stage of cascade	R 23

3rd stage for the cooling of treatment room

R 14

Width approx.:	2,400 mm
Depth approx.:	800 mm
Height approx.:	2,000 mm
Weight approx.:	1,360 kg
Ambient conditions:	
Min./max. ambient temperature:	5/25 °C
Sound pressure level at a distance of 3 meters, approx.:	72 dB (A)
Blow-off line cold:	1"/ min.PN 30bar (28 mm CU pipe)
Blow-off line water:	¾"/ min.PN 10bar (fabric hose/CU pipe)

Adequate ventilation and, if necessary, air conditioning of the machine room must be provided.

2.2 Electrical control cabinet, machine

Back-up fuse max.:	63 A
Maximum operating current at operating point*:	47.0 A
Total connected load at operating point*:	23 Kw
Electrical connection value:	3Ph/N/PE 400 V 50 Hz 63 A
Electrical system type:	TN – Network / TNC Network**
RCD installed:	Type A / 30 mA
Control unit:	Siemens PLC S7-1200
Cabinet:	Rittal
System Test:	accordance to DIN EN 60204-1
Relays / contactors:	Weidmüller / Siemens
Circuit breaker:	ABB
Main switch:	ABB
System test:	according to DIN-EN 60204-11
Width approx.:	800 mm
Depth approx.:	400 mm
Height approx.:	2,000 mm
Weight approx.:	167 kg

*The operating point is an empirically determined value, which can vary due to various external effects.

**If a residual current circuit breaker is used, it must be designed to be sensitive to universal current.

2.3 Electrical control cabinet, cell

The electrical control cabinet for the cell is powered via the electrical control cabinet for the machine. No power supply by the client necessary.

Cable routes from the machine electrical control cabinet to the cell electrical control cabinet are required.

Control unit:	Siemens PLC S7-1200
Cabinet:	Rittal
Relays / contactors:	Weidmüller / Siemens
Circuit breaker:	ABB
Main switch:	ABB
System test:	according to DIN-EN 60204-1

Width approx.:	1,200 mm
Depth approx.:	400 mm
Height approx.:	2,000 mm
Weight approx.:	200 - 243 kg

2.4 Special evaporator in the anteroom

- Evaporator designed specifically for wall or ceiling mounting
- Evaporator block soldered with silver solder
- Evaporator core tubes internally finned
- Fan specifically designed for low temperature applications
- Electric defrost incl. safety shut-down

2.5 Special evaporator in the treatment room

- Evaporator specifically designed for low temperature applications and wall mounting
- Evaporator block soldered with silver solder
- Evaporator core tubes internally finned
- Fan specifically designed for low temperature applications
- Electric defrost incl. safety shut-down
- 3 pcs external fan motors in special design with internal fan blades

2.6 Cell combination TYPE ISO 150/160

Insulation:	150 mm / 160 mm rigid polyurethane foam Foamed free of CFCs/HCFCs
Surface:	Sheet metal, galvanized on the inside and the outside, plastic coated
Pressure compensation:	4 pcs of pressure compensation valves, of which 2 with increased heating
Floor:	0,8 mm stainless steel floor manufactured as a bottom tray. Slip resistance class 11 load capacity 20kN/m ²
Lighting:	given in every room Außenabmessungen:
Exterior dimensions:	
Width:	4,200 mm
Depth:	2,400 mm
Height	2,550 mm
Weight approx.:	1,460 kg
OR:	
Exterior dimensions:	
Width:	4,200 mm
Depth:	2,100 mm
Height	2,800 mm
Weight approx.:	1,650 kg

Other external dimensions are available on request at extra charge.

At least 50 cm of space behind the cell is required for setting up the cell. If there is not enough space behind the cell, at least 50 cm is required above the cell.

The colour is similar to RAL 9010 on the inside and RAL 9006 on the outside, other RAL colours are possible, on request.

Environmental specifications:

Temperature range:	18-23 °C
Rel. air humidity, max.:	50 %

All joint seals on the cell, windows, doors and elements serve as a maintenance joint.

During the production of the cell elements it may occur outgassing in the foam, which can occur bubbles in the single cell elements. We exclude any warranty, because this is only an optical and not a technical defect.

2.6.1 All-glass entrance door

Type:	DIN left or right
Insulating glass door with frame	
Door dimensions approx.	926 x 2,045 mm (W x H)
Internal door dimensions approx.	869 x 1,991 mm (W x H)

- Circumferential magnetic seal
- Natural anodised aluminium frame with heating
- Continuous handle rail
- Heated triple glazing

2.6.2 Intermediate door

Type:	DIN left or right
Internal door dimensions approx.:	900 x 1900 mm (W x H)

- Intermediate door, single-leaf
- Temperature range -60 °C/-110 °C
- Foamed CFC-free, galvanized sheet steel similar to RAL 9010
- Designer handles for opening and closing, powder-coated, similar to RAL 9010
- Circumferential EPDM gasket sealing the frame
- Stainless steel door frame with special frame heating to prevent accumulation of ice
- Maintenance-free hinges and door clamps with aluminium rollers

2.6.3 Deluxe window to the treatment room

xtra-large window for the treatment room . Good visual contact between the operator and the user. Breaking up of the space, also suitable for users with a tendency for claustrophobia.

Dimensions: 500 x 1700 mm

2.7 Service router with interface for remote monitoring

- Internet connection via WAN port
- PLC monitoring with an alert system (via e-mail)
- VPN switch for simplified on-site administration
- German brand product

Network:

VPN	OpenVPN server, OpenVPN client
WAN/LAN	2 x RJ 45 10/100 Mbit/s, full and half duplex
WLAN	802.11 b/g/n, access point, WAN client

2.8 Camera in the treatment room

- Special camera incl. housing and window heater
- Special heat-resistant protective housing IP54
- PoE camera

2.9 Music playback in the treatment room

- 2 loudspeakers in treatment room

2.10 Communication system from control unit to treatment room

- IP-based
- Advanced audio technology, automatic gain and volume control, and background noise suppression

2.11 Operating unit and additional operating options

Display:	large 22 inch touch screen
Software:	Microsoft Windows 10

The following functions can be operated via the touch screen:

- Temperature monitoring

- Countdown timer
- Switching the light on and off
- Setting the timer


Remote VPN access can be booked optionally.

Example picture: Designer control panel made of carbon by capristo.de

Alternatively: stainless steel panel, chrome-plated in 2 different color variations:
Black or bronze chrome plated

Suitable for installation as wall-in mounting, as well as wall-on mounting.

For care instructions, please refer to point 6.

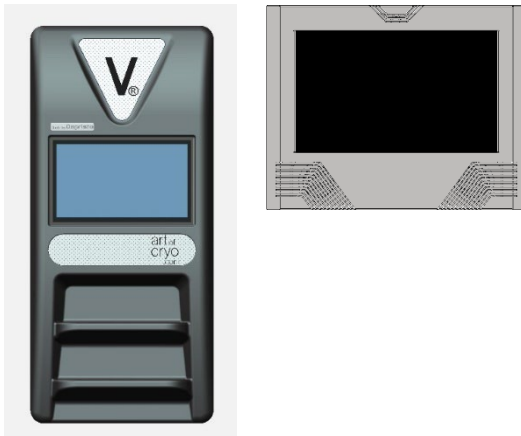2.12 Handrail in the treatment room

- Made of special hemlock wood

3. Optional components

- Water-cooled condenser ⚠ **no city-water (drinking-water)**
- Air-cooled condenser

3.1 Water-cooled condenser

⚠ **ATTENTION: no city-water (drinking-water)**

Intake/discharge connections:	1" R-AG
Cold water volume flow	3.5 m³/h
Water intake min./max.	5/20 °C
Water pressure min./max.	1.5/5 bar

Option: connection with town water, cold water set or with glycol (on request)


3.2 Air-cooled condenser

⚠ ATTENTION: no city-water (drinking-water)

- High-efficiency compact vertical condenser
- Housing: galvanized steel, similar to RAL 7035
- EC-fan which meets the efficiency requirements of Directive 2009/125/EC

Width approx.:	1,600 mm
Depth approx.:	800 mm
Height approx.:	1,400 mm
Weight approx.:	175 kg

Max. ambient temperature:	34 °C
Sound pressure level at a distance of 10 meters, approx.:	48 dB (A)

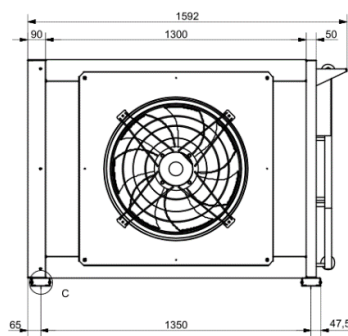4. Optional components

- Hybrid cooling (switchable air and water cooling)
- Emergency door
- Wood panelling incl. handrail in the treatment room
- Chiller
- Special dimensions V2 cell
- Big condenser
- Camera in the anteroom

- Music playback in the anteroom
- Heat recovery
- Division of the cooling unit
- Frost heave protection

4.1 Hybrid cooling

- Manually switchable air and water cooling on request for extra price

4.2 Emergency Door

Type:	DIN left or right
Door dimensions approx.	795 x 2,077 mm (W x H)
Internal door dimensions approx.:	869 x 1,992 mm (W x H)

- Temperature range -110 °C/+25 °C
- Foamed CFC-free, galvanized sheet steel similar to RAL 9010
- Designer handles for opening and closing, powder-coated, similar to RAL 9010
- Circumferential EPDM gasket
- Stainless steel door frame with special frame heating to prevent accumulation of ice
- Maintenance-free hinges and door clamps

4.3 Wood panelling incl. handrail in the treatment room

- Special hemlock wood panelling
- Cold-resistant
- Provides a pleasant ambience

4.4 Chiller

- Customised according to client requirements

4.5 Special dimensions V2 cell

- there are various options here, please contact us to find out more

4.6 Big condenser

For high ambient temperature up to:

40 °C

Width approx.:

2,950 mm

Depth approx.:

770 mm

Height approx.:


1,450 mm

Sound pressure level at a distance of 10 meters, approx.:

56 dB (A)

Weight approx.:

278 kg

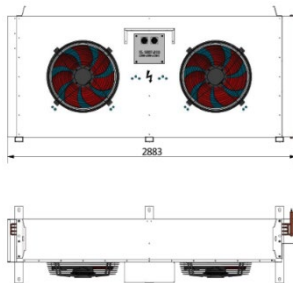4.7 Camera in the anteroom

- Special camera incl. housing and window heater
- Special heat-resistant protective housing IP54
- PoE camera

4.8 Music playback in the anteroom

- Loudspeaker: 2-way Hi-Fi coaxial loudspeaker
- Amplifier: 2x 15 / 10 W / 4 Ohm / 8 Ohm or 1x 30 W / 8 Ohm; power supply 12-24 V incl. convection cooling
- Coupling through Bluetooth receiver

4.9 Heat recovery

The condensation capacity is 16 kW.

- Delta 5 K
- Connection nozzle 1 inch
- Pressure drop approx. 0.3 bar
- Flow rate 2.76 m³

The heat output depends on the desired water temperatures. High water temperatures lead to partial liquefaction in the heat recovery heat exchanger, the remaining liquefaction takes place in the downstream condenser.

Example with a condensation power of 16 kW

100% heat recovery at a water temperature of up to 38 °C:	approx. Q = 16.0 kW
approx. 60% heat recovery at a water temperature of up to 45 °C:	approx. Q = 13.0 kW
approx. 40% heat recovery at a water temperature of up to 50 °C:	approx. Q = 9.0 kW
approx. 20% heat recovery at a water temperature of up to 55 °C:	approx. Q = 4.0 kW

Information: When using chlorinated water in operation, an intermediate circuit must be installed to separate the system. A titanium heat exchanger can be used as an option (on request).

4.10 Division of the cooling unit

- The unit can be divided in width for better set-up, e.g. in narrow rooms (for extra price)

4.11 Frost heave protection

Our low-temperature rooms give off cold to their surroundings. If this takes place next to foundations that are partially or completely embedded in the ground, the ground picks up the cold. If the natural heat flow from the environment and the soil is not sufficient, the moisture in the soil will freeze, expand and heave the soil up. Depending on the soil conditions, this can cause considerable damage to the building. The use of frost heave protection must be verified before installation.

5. On-site services

In order to ensure optimal introduction, installation and commissioning of V2, please meet the following prerequisites:

5.1 Connections

1. Current:

Back-up fuse max.:	63 A
Operating voltage:	3Ph/N/PE 400 V 50 Hz
Electrical system type:	TN - Network / TNC Network

2. Network / Internet:

The service router sets up an outbound TCP connection to the address.

<https://vpn.kaeltecloud.de> (IPv4: 128.140.2.254, alt. IPv6: 2a01:4f8:c2c:3e7f::1) In order to establish a secure connection to the server, the following requirements must be met:

- a DNS server must be accessible in your network
- the above IP address (IPv4: 128.140.2.254) must be enabled
- the following TCP-Port must be enabled for the VPN tunnel: TCP-Port 443
- the following TCP-Port must be enabled for the Cloud-function: TCP-Port 8883
- UDP port 123 should be open for synchronisation

The system is fitted with an RCD Type A / 30 mA.

5.2 Set-up

- Unloading and bringing in of the components
- Container or similar for the disposal of the packaging material
- Free routes for bringing in
- If necessary floors should be protected (no warranty for damage)
- Fire safety measures
- Masonry, chiselling, earthworks, roofing and coating
- Lifting equipment such as cranes, scaffoldings, mobile working platforms, forklift trucks etc.

- Installation on level floor with flowing screed/ tiles, 2 cm tolerance over a total length of 6 m
- Breakthrough as hole for condenser/glycol incl. drain line: 200 mm diameter (for all variants).

Breakthrough between machine and cell:

single Vaultz: 200 x 500 mm

vario Vaultz: 300 x 700 mm

All breakthroughs must be sealed by the customer after our installation according to local fire standards.

In order to guarantee problem-free function of V2, the following points must be taken into account:

5.3 Ambient conditions

1. Installation room cell:

- Around the cell at least 5-fold air exchange of the room volume according to DIN EN 16798 part 3.
- Temperature range: 18-23 °C
- Rel. humidity: max. 50%
- If a dehumidifier is used, please ensure automatic water drainage.
- It is possible that moisture (e.g. condensation or water due to snow formation) may escape from the cell. The manufacturer accepts no liability for any damage to the substrate and/or parts such as walls or drywall components. Please take appropriate precautions here.

2. Machine room:

Sufficient intake and exhaust ventilation and, possibly, air conditioning of the machine room, to make sure that an ambient temperature of min./max.: 5 / 25 °C is maintained.

Please note that

a) the heat emission of the cooling unit and control cabinets is max. 8 kW

b) There are other heat sources and a buffer.

c) We expressly point out that the machine can transmit disturbing vibrations and noises to the building during operation. If necessary, precautions against the transmission of vibrations and noises to the building must be organized by the customer. This is not the responsibility of the manufacturer.

Notes when using a container as a machine room:

- Up to max. 22°C outside temperature, the machine room can be ventilated and de-aerated via the fans and the ventilation grids

(Note: the ventilation grilles must be regularly checked for dirt and cleaned if necessary)

- From 22°C outside temperature, fans must be switched off via the thermostat and the ventilation grilles must be closed.

- Air-conditioning must be switched on from 22°C

(Note: the air conditioning must be checked regularly for proper functioning)

6. General provisions

- Design and documentation of the machine according to the L&R standard and CE certification
- The limit values specified in EN 50160 must be observed for supply voltage

Note on water quality:

Various processes can take place in water-bearing systems, which can lead to corrosion, mineral deposits and also biological problems. We reject any claims for malfunctions or damage to systems or system components caused by insufficient water quality. In order to achieve proper and trouble-free operation of water-bearing systems, the water quality must therefore be checked and, depending on system-relevant data, water treatment and water treatment must be carried out.

On request, we can provide you with a reference value table for the recommended water quality.

Care instructions for the stainless steel operating unit:

Please note that aggressive cleaners must not be used. We recommend using a chrome cleaner with a care function. We also recommend using a microfiber cloth.

Technical changes:

Our systems, machines and components are state-of-the-art. We reserve the right to introduce technical changes that do not impair the performance and functions of the system.

Cable lengths

Machine to cell	= 10 m
Machine to condenser	= 10 m

General note:

The shorter the pipes between the machine and the cell, the more energy-efficient the machine.

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EU Type examination (module B) - production type - according to Directive 2014/68/EU

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CERTIFICATE

International Certification Management GmbH

L&R Kältetechnik GmbH & Co. KG
Innovative Kältetechnik
D-59846 Sundern Hachen, Hachener Str. 90a

Standard: **DIN EN ISO 9001:2015**

Scope:
Development, installations, distribution of refrigeration systems, switch and control cabinets, refrigeration and cryogenic systems for process engineering and manufacturing operations in industry

With this certificate it is confirmed, that the L&R Kältetechnik GmbH & Co. KG Innovative Kältetechnik applies a management system which meets the requirements of DIN EN ISO 9001:2015

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Issue date
02.12.2025

Valid until
01.12.2028

Translation of the German original

International Certification Management GmbH
Straubinger Straße 19
94363 Oberschneiding
Germany

info@ic-management.com
www.ic-management.com

International Certification Management

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by L&R Kältetechnik GmbH & Co. KG